




# Accessories

## Telescopic arm and swivelling divider

Dividing curtain specially designed to preserve privacy between patients, being able to monitor them by a standing up person.

- The support is extended by pulling the control located on one of the edges and is gathered again by pushing in the opposite direction, being possible also to arrange it in parallel to the wall to clear room when it is not in use.
- Rail version, wall version and Ø38,0 mm tube version available.

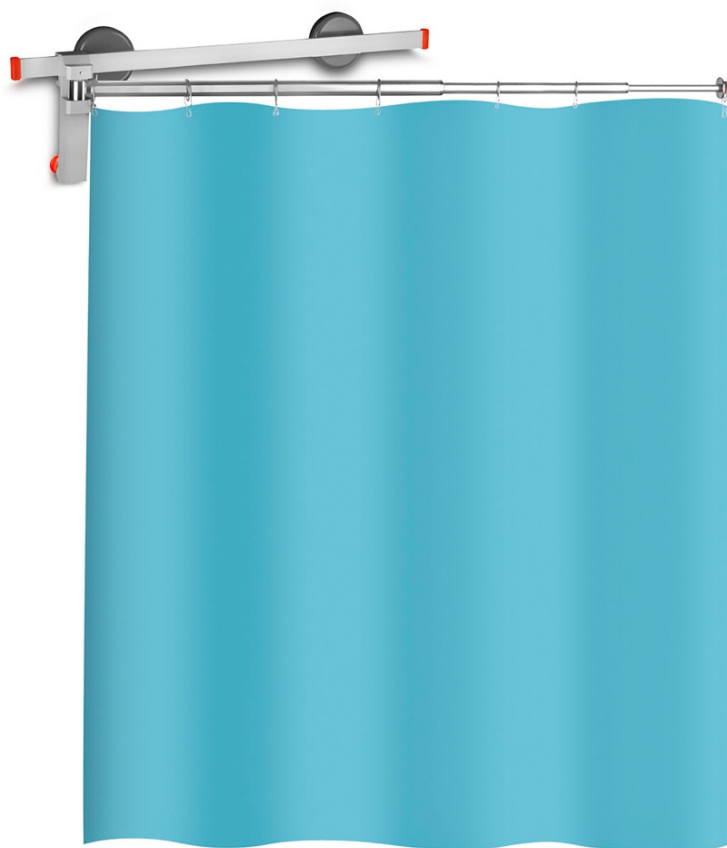

### TECHNICAL SPECIFICATIONS

- Manufactured in anodized aluminium, chromium-plated brass and stainless steel.
- Length gathered / extended: 590 / 1510 mm.
- Curtain area: 1500 x 1000 mm.
- Delivered in a package of 5 disposable curtains.